




CARLO Comfort GP L

Distributed by
PRISM MEDICAL



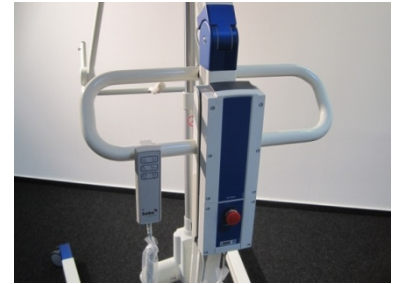

Ergonomic handles



Lockable carry bar



Hand control

Aluminium lift CARLO Comfort GP L

... the multi-purpose shower and transfer lift

The Carlo GP L utilizes a gas strut making it easy to sit your client upright or recline for comfort. With the CARLO Comfort GP L, transferring to a bed or chair is made easy. Features an aluminum powder-coated chassis with electrically controlled spreading legs. Maneuvering this lift is made easy with 4 double closed system castors. The 2 rear castors are lockable. The CARLO Comfort GP L is a unique floor based lift that easily integrates with any clip style sling. The unique design maximizes space for the client allowing for lifts from bed, chair or floor. It is available in a variety of configurations including: digital scales, dual controls, on-board charging or a dual battery pack system. The Carlo can be configured to meet the needs of any facility.

Technical data:

Material: Aluminium - anodized/powder coated
Base: (length/height) 129 cm / 12 cm
Leg width inside: (closed/open) 58 cm / 104 cm
Leg width outside: (closed/open) 70 cm / 114 cm
Total width of the lift (base legs closed): 70 cm
Lifting arm range: (low/high) 68 cm / 208 cm
Maximum patient weight: 181.5 kg (400 lbs)
Power supply: 24 volt Linak battery
Lifting motor: 24 volt Linak with manual emergency lowering

CARLO Comfort GP L Part # 278850
(on board charging and dual controls)

CARLO Comfort GP L Part # 278855
(includes; 2 x batteries, wall mount charger and dual controls)

Accessories:

CARLO Comfort GP L Scale Part # 278801

Transfer slings, bathing slings, Toilet slings are available in all sizes

Am Rübenmorgen 3
D-35582 Wetzlar-Dutenhofen
Fon +49 (0) 641/92220-0
Fax +49 (0) 641/92220-20
Email: info@beka-hospitec.de
www.beka-hospitec.de



Distributed by Prism Medical Ltd.

480 University Avenue, Suite 100, Toronto, ON Canada M5G 1V2

T 1.800.265.0677 • www.prismmedical.ca

PRISM MEDICAL